



CAMERAS & DIGITAL Solutions

Cameras for Microscopy & Digital Scanner


CAMERAS & DIGITAL Solutions

4

USB CAMERAS - B Series	page 375
USB CAMERAS & TABLETS PC - TB Series	page 378
USB CAMERAS - PRO Series	page 380
PRO Series Cooled	page 384
HDMI CAMERAS - H Series	page 386
HDMI CAMERAS & SCREEN - H Series	page 388
HDMI AUTO-FOCUS CAMERAS - C-HA	page 390
WiFi & USB CAMERA - 4083.WiFi	page 392
EYEPIECE CAMERAS - Educam & VC Series	page 394
OPTICAL ADAPTERS - Adapters for OPTIKA Cameras and Microscopes	page 396
OPTIKA SOFTWARE - Microscopy Analysis Software Suite	page 397
OPTISCAN - Digital Slide Scanner	page 405



Icons

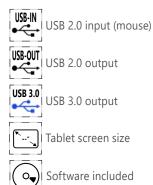



4

Microscope Cameras




MICROSCOPE CAMERAS

USB CAMERAS - B Series Series Output Description: Descriptindescriptindescription:

C-B1 / C-B3 / C-B5 / C-B10



USB user-frienldy cameras for **general purposes** Superb results and vivid details from standard to high resolution

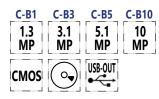

USB 2.0 C-Mount and Eyepiece Microscope Cameras For General Purposes


Cameras have become indispensable nowadays and OPTIKA is offering a line of remarkable solutions for digital imaging. OPTIKA B Series represents a cost-effective solution equipped with the latest technology sensors with more vivid colors and great contrast for stunning images. This series features Aptina CMOS sensor with excellent color reproduction, significantly high frame rates and several resolutions available to match any customer need. Thanks to the convenience and simplicity, being extremely intuitive to install and operate, the OPTIKA B Series is recommended for educational and routinary microscopes, also as eyepiece cameras (no need for additional adapters/rings in case of monocular and binocular microscopes). Software included (for Windows): Optika Vision Lite, Optika LiteView and Optika ProView. Optika LiteView is also compatible with IOS & Linux.

USB CAMERAS - B Series - Specifications

	С-В1	С-ВЗ
Digital camera resolution	1.3 MP (1280 x 1024)	3.1 MP (2048 x 1536)
Signal output	USB 2.0	USB 2.0
Sensor Size	1/3″	1/2"
Sensor technology	CMOS	CMOS
Sensor type	Aptina CMOS	Aptina CMOS
Image format	5/4	4/3
Pixel size	3.6 x 3.6 µm	3.2 x 3.2 μm
Frame rate full resolution	15 fps (1280 x 1024)	12 fps (2048 x 1536)
Frame rate other resolutions	50 fps (320 x 256)	32 fps (1024 x 768); 45 fps (680 x 510)
Sensitivity	1 V/lux-second	1 V/lux-second
Signal / noise ratio	44 dB	43 dB
Dynamic range	74 dB	61 dB
ADC conversion	8 Bit	8 Bit
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit	1 Bit ; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.14 msec - 2 sec	0.244 msec - 2 sec
Binning	1x1; 2x2; 4x4	1x1; 2x2; 3x3
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
Camera power	PC USB	PC USB
C-mount	YES	YES

Accessories included: C-mount projection lens, 30 mm / 30.5 mm ring adapters, calibration slide, 1.8 m USB cable.

USB CAMERAS - B Series

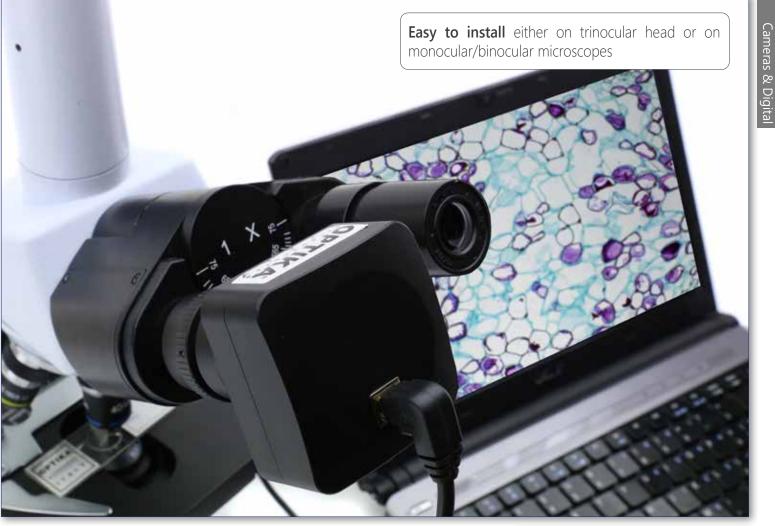

С-В5	C-B10
5.1 MP (2592 x 1944)	10 MP (3584 x 2748)
USB 2.0	USB 2.0
1/2.5″	1/2.3″
CMOS	CMOS
Aptina CMOS	Aptina CMOS
4/3	4/3
2.2 x 2.2 μm	1.67 x 1.67 μm
7 fps (2592 x 1944)	3.3 fps (3584 x 2748)
27 fps (1280 x 960); 90fps (640 x 480)	11 fps (1792 x 1374); 38 fps (896 x 684)
0.53 V/lux-second	0.31 V/lux-second
40.5 dB	34 dB
66.5 dB	65.2 dB
8 Bit	8 Bit
1 Bit; 4 Bit; 8 Bit; 24 Bit	1 Bit; 4 Bit; 8 Bit; 24 Bit
0.294 msec - 2 sec	0.4 msec - 2 sec
1x1; 2x2; 4x4	1x1; 2x2; 4x4
380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
PC USB	PC USB
YES	YES



4 **USB CAMERAS & TABLETS PC** - TB Series





Unique Features

- » Simultaneous camera & power connection
- » Equipped with the latest Windows OS & Intel processor
- » Easily detachable, can be used as a laptop (keyboard included)

USB 2.0 C-Mount Microscope Cameras with Windows Tablet PC For **Unparalleled Comfort & Performance**

Exclusive tablet PC, powerful and versatile for a great user experience. Always one step forward to ensure the latest technology! • A 2-in-1 solution that you can use like a PC, being Windows-based

- Powerful Intel processor ensuring top performance and speed
- High-resolution, vivid color graphic display
- Large touch screen of 10.1" with fast, responsive and smooth control
- Attached camera available in 3.1 MP (TB-3W) or 5.1 MP (TB-5W) resolution
- Includes the user-friendly and intuitive Optika ProView software
- C-mount connection for trinocular microscopes only

TB-3W





Windows tablet PC with large 10.1" LCD touch screen, combined with a 3.1 MP camera to create the most advanced solution for digital microscopy.



Windows tablet PC with large 10.1" LCD touch screen, combined with a 5.1 MP camera to create the most advanced solution for digital microscopy.

USB CAMERAS & TABLETS PC - TB Series



Perform linear measurement on your image with **Optika ProView** just by drawing a line!

Cameras & Digital



TABLET TECHNICAL SPECIFICATIONS		
Operating system	Windows 10 (64Bit)	
CPU	Intel® Atom™ Z8350, Quad core	
CPU speed	1.44 GHz	
Graphic card	Intel® HD Graphics 400	
RAM	Ram 4 GB LPDDR3	
Display size	LED 10.8" IPS Multi Touch Screen	
Display resolution	1920x1200	
Storage	Hdd 64 GB	
Network	Wireless - Bluetooth 4.0	
Input ports	Micro USB - USB - Micro SD card reader	
Output ports	Microphone - Headphone - Micro HDMI	
Battery Type	Lithium-ion	
Battery capacity	6500 mAh	
Power consumption	15W	
Power supply	5V 3A	
Dimensions (mm)	261 X 167 X 9 mm	
Weight (Kg)	0.53 kg	
Language	Multilanguage	

CAMERA TECHNICAL SPECIFICATIONS	TB-3W	TB-5W
Digital camera resolution	3.1 MP (2048 x 1536)	5.1 MP (2592 x 1944)
Signal output	USB 2.0	USB 2.0
Sensor Size	1/2″	1/2.5″
Sensor technology	CMOS	CMOS
Sensor type	Aptina CMOS	Aptina CMOS
Image format	4/3	4/3
Pixel size	3.2 x 3.2 μm	2.2 x 2.2 μm
Frame rate full resolution	12 fps (2048 x 1536)	7 fps (2592 x 1944)
Frame rate other resolutions	32 fps (1024 x 768); 45 fps (680 x 510)	27 fps (1280 x 960); 90fps (640 x 480)
Sensitivity	1 V/lux-second	0.53 V/lux-second
Signal / noise ratio	43 dB	40.5 dB
Dynamic range	61 dB	66.5 dB
ADC conversion	8 Bit	8 Bit
Color Depth	1 Bit ; 4 Bit; 8 Bit; 24 Bit	1 Bit; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.244 msec - 2 sec	0.294 msec - 2 sec
Binning	1x1; 2x2; 3x3	1x1; 2x2; 4x4
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
Camera power	PC USB	PC USB
C-mount	YES	YES

Accessories included: C-mount projection lens, calibration slide, 1.8 m USB cable , keyboard with touchpad and touch pen.

⁴ USB CAMERAS - PRO Series

C-P3 / C-P6 / C-P8



Professional cameras with high sensitivity, low noise and impressive capabilities in special applications (on different observation methods)

C-P3	C-P6	C-P8
3.1 MP	6.3 MP	8.3 MP
CMOS	•	USB-OUT

High-Performance USB 3.0 C-Mount Microscope Cameras For General Purposes With High FPS

Do you require a high-end camera with an especially high resolution, generous dynamic range, rapid read-out rate and a USB3.0 port?

If your answer is yes, then the PRO series is your choice. Its compact and elegantly designed housing conceals the very latest in camera technology. Your images will be of the highest quality and rich in contrast and detail. OPTIKA PRO Series includes a wide range, to virtually fulfill each application demand.

Top-class SONY sensors, worldwide recognized, ensure you to capture your specimen in beautiful true-to-life color, delivering incredibly accurate colors just as you see them. Software included (for Windows): Optika Vision Lite, Optika LiteView and Optika ProView. Optika LiteView is also compatible with IOS & Linux.

USB CAMERAS - PRO Series - Specifications

	С-РЗ	C-P6	С-Р8
Digital camera resolution	3.1 MP (2048 x 1536)	6.3 MP (3072 x 2048)	8.3 MP (3840 x 2160)
Signal output	USB 3.0	USB 3.0	USB 3.0
Sensor Size	1/2.8″	1/1.8″	1/2.5″
Sensor technology	CMOS	CMOS	CMOS
Sensor type	SONY EXMOR CMOS	SONY EXMOR CMOS	Sony exmor cmos
Image format	4/3	3/2	16/9
Pixel size	2.5 x 2.5 µm	2.4 x 2.4 µm	1.62 x 1.62 μm
Frame rate full resolution	50 fps (2048 x 1536)	30 fps (3072 x 2048)	32 fps (3840 x 2160)
Frame rate other resolutions	50 fps (1920 x 1080)	38 fps (1536 x 1024)	65 fps (1920 x 1080)
Sensitivity	600 mV at 1/30sec	425 mV at 1/30sec	236 mV at 1/30sec
Dark signal	0.15 mV at 1/30sec	0.15 mV at 1/30sec	0.1 mV at 1/30sec
ADC Conversion	8 Bit - 12 Bit	8 Bit - 12 Bit	8 Bit - 12 Bit
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit	1 Bit; 4 Bit; 8 Bit; 24 Bit	1 Bit; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.1 msec - 15 sec	0.1 msec - 15 sec	0.1 msec - 15 sec
Binning	1x1	1x1; 2x2	1x1; 2x2
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
Camera power	PC USB	PC USB	PC USB
C-mount	YES	YES	YES

Accessories included: calibration slide, 1.8 m USB cable

USB CAMERAS - PRO Series

4083.F33



USB 3.0 C-Mount Microscope Camera With High FPS

Excellent color reproduction, with high frame rates, the OPTIKA 4083.F33 represents a valuable choice in routinary microscopy. Software included (for Windows): Optika Vision Lite



4083.F33
3 MP
USB 3.0
1/2.7"
CMOS
ISP Integrated Chip
16/9
3.0x3.0 μm
30 fps (1920x1080)
7 fps (2304x1936)
1.275 V/lux-sec (550 nm)
na
8 Bit
8 Bit
330 msec
1x1; 2x2
700 nM
PC USB
YES
Accessories included: calibration dide 1.9 m LISP cable

Accessories included: calibration slide, 1.8 m USB cable.

Achieve professional results and outstanding images with PRO Series: incredibly sharp and high contrast images also when combined to a stereomicroscope

ISB CAMERAS - PRO Series

C-P6FL



High-performance CCD camera for special applications. Recommended for **professional microscopes**, routinary and research level (on different observation methods)



High-Performance USB 3.0 C-Mount Microscope Cameras For Special Applications

High grade SONY EXView HAD CCD sensor camera with 6.0 MP resolution recommended for special applications thanks to the impressive capabilities in the most challenging working conditions and featuring high sensitivity, low noise and USB 3.0 connection.

Software included (for Windows): Optika Vision Lite, Optika LiteView and Optika ProView. Optika LiteView is also compatible with IOS & Linux.



USB CAMERAS - PRO Series - Specifications


	C-P6FL	
Digital camera resolution	6 MP (2748 x 2200)	
Signal output	USB 3.0	
Sensor Size	1″	
Sensor technology	CCD	
Sensor type	SONY ExView HAD	
Image format	5/4	
Pixel size	4.54 x 4.54 μm	
Frame rate full resolution	7.5 fps (2748 x 2200)	
Frame rate other resolutions	14 fps (2748 x 1092)	
Sensitivity	1000 mV at 1/30sec	
Dark signal	8 mV at 1/30sec	
ADC Conversion	8 Bit - 14 Bit	
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit	
Exposure Time	0.06 msec - 1000 sec	
Binning	1x1	
IR filter	380-650 nm (IR-cut filter)	
Camera power	PC USB	
C-mount	YES	

Accessories included: calibration slide, 1.8 m USB cable.

Several focusable C-mount adapters available



USB CAMERAS - PRO Series

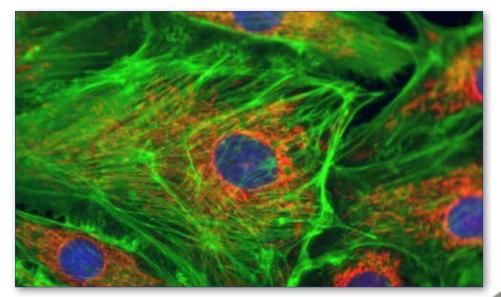

L10X/22

CCD sensor ensures significant and **impressive quality** on special techniques, **including fluorescence**.

Non-cooled high-performance CCD camera. Thanks to its very high sensitivity, it's perfect for fluorescence microscopy.

USB CAMERAS - PRO Series Cooled

C-P20CC / C-P20CM





Ultra High-Performance USB 3.0 C-Mount Cooled Microscope Cameras For Special Applications

For special applications, and requiring different observation methods, high sensor size cameras are preferred for the high sensitivity delivering high signal-to-noise ratio (low noise is achieved via cooling) and a large dynamic range. OPTIKA Pro Cooled cameras provide an excellent sensitivity. Peltier-cooled (cooling to 45°C below ambient), with scientific grade CMOS sensor, they ensure great performance in low light conditions, ultra-long exposure time and an impressive reliable color fidelity.

Software included (for Windows): Optika Vision Lite, Optika LiteView and Optika ProView. Optika LiteView is also compatible with IOS & Linux.



USB CAMERAS - PRO Series Cooled - Specifications


	C-P20CC (Color)	C-P20CM (Monochrome)
Digital camera resolution	20 MP (5440 x 3648)	20 MP (5440 x 3648)
Signal output	USB 3.0	USB 3.0
Sensor Size	1″	1″
Sensor technology	CMOS	CMOS
Sensor type	SONY EXMOR	SONY EXMOR
Image format	3/2	3/2
Pixel size	2.4 x 2.4 µm	2.4 x 2.4 µm
Frame rate full resolution	5 fps (5440 x 3648)	17.8 fps (5440 x 3648)
Frame rate other resolutions	10 fps (4096 x 2160); 15 fps (2736 x 1824); 30 fps (1824 x 1216)	41 fps (4096 x 2160); 51 fps (2736 x 1824); 64 fps (1824 x 1216)
Sensitivity	426mV at 1/30sec	388mV at 1/30sec
Dark Signal	0.21mV at 1/30sec	0.21mV at 1/30sec
Cooling System	Internal two-stage TE cooling system -45°C	Internal two-stage TE cooling system -45°C
ADC conversion	8 Bit - 14 Bit	14 Bit
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit	1 Bit; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.1 msec - 3600 sec	0.1 msec - 3600 sec
Binning	1x1; 2x2; 3x3	1x1; 2x2; 3x3
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
Camera power	PC USB	PC USB
C-mount	YES	YES

USB CAMERAS - PRO Series Cooled





Durable and safe reinforced plastic case to protect and store your valuable camera

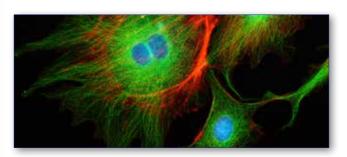

Achieve unparalleled results and superb images



IDMI CAMERAS - H Series



720p HD high-speed microscope camera with Aptina CMOS sensor: intuitive, entry-level

C-HE




C-Mount & Eyepiece HDMI Microscope Camera For Easy Operation On Screens, No Compromises In Quality

- Aptina CMOS sensor with excellent color reproduction for routine applications
- C-mount connection, eyepiece adapter available on request
- No installation of software required when used in HDMI mode
- HDMI (2 MP, 1280 x 720) camera for screen connection
- Extremely reliable color fidelity
- Lightweight cameras, to be used even on the smallest and lightest microscopes
- SD card enables image and video capturing
- Built-in function buttons for HDMI camera control



C-HP

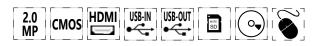

C-Mount & Eyepiece High-Performance HDMI / USB Microscope Camera For Professional Use On Screens & PC, With Maximum Flexibility

- High grade SONY CMOS sensor for reliable colors
- C-mount connection, eyepiece adapter available on request
- No installation of software required when used in HDMI mode
- HDMI FULL HD (2 MP, 1920 x 1080) camera for screen connection and USB camera for PC use
- Reliable color fidelity
- Highly reccomended for wide range of applications
- SD card enables image and video capturing
- Software included (for Windows): Optika ProView





Professional 1080p FULL HD high-speed microscope camera with measurement functions for **premium performance**

4

HDMI CAMERAS - H Series



HDMI CAMERAS - H Series - Specifications

	C-HE	С-НР
Video resolution (USB output)	-	1920 x 1080 pixel
Video resolution (HDMI output)	HD 720p	Full HD 1080p
Digital camera resolution	2 MP (1280 x 720)	2 MP (1920 x 1080)
Signal output	HDMI	HDMI
Sensor Size	1/2.8″	1/1.9″
Sensor technology	CMOS	CMOS
Sensor type	APTINA	SONY
Image format	16/9	16/9
Pixel size	2.8 x 2.8 µm	3.75 x 3.75 µm
Frame rate (HDMI)	30 fps (1280 x 720 HDMI); 30 fps (1920 x 1080 Capture)	60 fps (1920 x 1080 HDMI); 26fps (1920 x 1080 USB)
Sensitivity	510 mV at 1/30sec	1120 mV at 1/30sec
Dark Signal	0.15mV at 1/30sec	0.15mV at 1/30sec
Exposure Time	0.06 msec - 1900 msec	0.34 msec - 4 sec
Binning	1x1	1x1
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
Camera power	DC 12V 1A	DC 12V 1A
C-mount	YES	YES
White balance	Auto	Auto / Manual
Gain control	Manual	Auto / Manual
Exposure control	Auto / Manual	Auto / Manual

Accessories included: C-HE: 2 m HDMI cable, 16GB SD card, wireless mouse, power supply.

C-HP: Calibration slide, 1.8 m USB cable, 2 m HDMI cable, 16GB SD card, wireless mouse, power supply.

HDMI CAMERAS & SCREEN - H Series C-HESC





C-Mount HDMI Microscope Camera With 11.5" Full HD LCD Screen, No Compromises In Quality

- All-in-one: save space on the bench, with no need for external devices
- Feel free to tilt the screen and adjust the position to the preferred level, eliminating fatigue during observation
- Aptina CMOS sensor with excellent color reproduction for routine applications
- HDMI (2 MP, 1280 x 720) camera
- C-mount connection for trinocular microscopes only



C-HPSC



C-Mount High-Performance HDMI / USB Microscope Camera with 11.5" Full HD LCD Screen For Professional Use (also on PC), With Maximum Flexibility

- All-in-one: save space on the bench, with no need for external devices
- Feel free to tilt the screen and adjust the position to the preferred level, eliminating fatigue during observation
- High grade SONY CMOS sensor for reliable colors
- HDMI FULL HD (2 MP, 1920 x 1080) camera for screen connection and USB camera for PC use
- Camera control panel shows exposure, gain, white balance, color adjustment, sharpness and noise control
- On-screen monitor toolbar shows measuring, mirror, comparison, zoom, freeze, cross & browser when using mouse control
- Software included (for Windows): Optika ProView
- C-mount connection for trinocular microscopes only





HDMI CAMERAS & SCREEN - H Series



HDMI CAMERAS & SCREEN - H Series - Specifications

	C-HESC	С-НРЅС
Video resolution (USB output)	-	1920 x 1080 pixel
Video resolution (HDMI output)	HD 720p	Full HD 1080p
Digital camera resolution	2 MP (1280 x 720)	2 MP (1920 x 1080)
Signal output	HDMI	HDMI
Sensor Size	1/2.8″	1/1.9″
Sensor technology	CMOS	CMOS
Sensor type	APTINA	SONY
Image format	16/9	16/9
Pixel size	2.8 x 2.8 µm	3.75 x 3.75 μm
Frame rate (HDMI)	30 fps (1280 x 720 HDMI); 30 fps (1920 x 1080 Capture)	60 fps (1920 x 1080 HDMI); 26fps (1920 x 1080 USB)
Sensitivity	510 mV at 1/30sec	1120 mV at 1/30sec
Dark Signal	0.15mV at 1/30sec	0.15mV at 1/30sec
Exposure Time	0.06 msec - 1900 msec	0.34 msec - 4 sec
Binning	1x1	1x1
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
Camera power	DC 12V 1A	DC 12V 1A
C-mount	YES	YES
White balance	Auto	Auto / Manual
Gain control	Manual	Auto / Manual
Exposure control	Auto / Manual	Auto / Manual
MONITOR TECHNICAL SPECIFICATIO	ons	
Size	11,5″	

Size	II,5"
Power supply	12V / 2,5A
HDMI cable	50 cm

Accessories included:

C-HE: 2 m HDMI cable, 16GB SD card, wireless mouse, power supply.
 C-HP: Calibration slide, 1.8 m USB cable, 2 m HDMI cable, 16GB SD card, wireless mouse, power supply.

neras & Digi

ta

Momentation Auto-Focus Cameras - C-HA C-HA



Equipped with the **state-of-the-art autofocus** system to ensure precise focusing in any condition and in **real time**



C-Mount High-Performance Autofocusing HDMI Microscope Camera

A superb microscope camera with the state-of-the-art autofocusing system, ensuring precise and ultra-fast automatic focus adjustment in any condition and in real time. Ideal to compensate the lack of parfocality of the microscope without any user effort.

- Aptina CMOS sensor with excellent color reproduction for routine applications
- HDMI FULL HD (2 MP, 1920 x 1080) camera for screen connection
- Camera control panel shows exposure, white balance, color adjustment and sharpness when using mouse control
- C-mount connection for trinocular microscopes only



HDMI AUTO-FOCUS CAMERAS - C-HA - Specifications

	С-НА
Digital camera resolution	2 MP (1920 x 1080)
Signal output	HDMI
Sensor Size	1/2.8″
Sensor technology	CMOS
Sensor type	APTINA CMOS
Image format	16/9
Pixel size	2.9 x 2.9 μm
Frame rate (HDMI)	50 fps (1920x1080)
Sensitivity	510 mV at 1/30 sec
Dark Signal	0.15 mV at 1/30 sec
Exposure Time	Auto
Binning	1x1
IR filter	YES
Camera power	5V 2A
C-mount	YES
CS-Mount	YES

Accessories included: 2 m HDMI cable, 16GB SD card, wireless mouse, power supply.



HDMI AUTO-FOCUS CAMERAS - C-HA





Compensates any lack of parfocality of the objectives

WIFI & USB CAMERA - 4083.WIFI

4083.WiFi





WiFi & USB 2.0 C-Mount and Eyepiece Microscope Camera Ideal For Teaching & Discussion Groups

- WiFi & USB connection modes
- 5 MP resolution on USB mode and 2 MP resolution on WiFi mode
- Ideal for educational applications with direct browser / IP connection
- No router is required
- Connectable to PC, smartphones, tablets and any other device
- Approx. 250 users connectable with good speed (speed depends on users connected)
- Supported by any device (PC, tablet or smartphone) with any type of browser (Chrome recomended)
- Lightweight cameras, to be used even on the smallest and lightest microscopes
- Image and video capturing function when used in WiFi mode
- Connection for both trinocular C-mount port and eyepiece tube
- Software included (for Windows): Optika Vision Lite, on USB mode

WIFI & USB CAMERA - 4083.WIFI - Specifications

	OPTIKAM WIFI 4083.WiFi
PC camera resolution	5 MP
WiFi camera resolution	2 MP
Signal output	USB 2.0, WiFi
Sensor Size	1\2.5″
Sensor technology	CMOS
Image format	4\3
Full Image size	2592 x 1944
USB Frame rate Full resolution	3 fps (2592 x 1944)
USB Frame rate other resolutions	11 fps (640 x 480) 8 fps (1024 x 768)
WiFi Frame rate Low resolution	2 fps (640 x 480) with 10 users
WiFi Frame rate other resolutions	1 fps (1024 x 768) with 10 users 1 fps/3sec: (1600 x 1200) with 10 users
WiFi Speed	802.11n 150Mbps
WiFi Antenna	2.4GHz 5DB
Max Exposure time	Automatic
External camera power	5 V 2000mA
White Balance	Auto / Man
Gain Control	Auto / Man
Back light control	Auto / Man
Exposure control	Auto / Man
C-Mount connection	YES
CS-Mount connection	YES, ready

Accessories included: C-mount projection lens, 30 mm / 30.5 mm ring adapters, calibration slide, 1.8 m USB cable.

USB & WiFi camera in once. No router is required



WIFI & USB CAMERA - 4083.WIFI



4 **EYEPIECE CAMERAS** - Educam & VC Series

Educam / VC-05



EDUCAM - Multimedia cameras to meet various requirements in the educational field

- Direct connection to TV screen and monitor
 - Versatile and flexible, yet sturdy and stable at the same time: can be used as overhead projector, for the projection of drawings, as a camera for teleconferences, assemblies, meetings or as a camera for filming
- Up to 90x magnifying power for any specimen and object
- 8mm objective lens enables focus from 0,76 cm, up to an infinite distance
- Extremely sensitive microphone to record voices/sounds (Multimedia models only)
 - Includes the user-friendly and intuitive OPTIKA Vision Lite software

All models are equipped with two adapters for video-microscopy (for biological and stereo microscopes)

VC-05 - Simple Eyepiece Microscope Camera With CCD Sensor, 420 TV Lines (PAL)







4083 - 4083.1 - 4083.2 - 4083.3

USB - 4083.4

MIC-4083.5

	MULTIMEDIA / 4083	MULTIMEDIA PRO / 4083.1	STUDENT / 4083.2	STUDENT PRO / 4083.3
Digital camera resolution	NO	NO	NO	NO
Analog camera resolution	PAL 582 x 420	PAL 582 x 420	PAL 582 x 420	PAL 582 x 420
Signal output	PAL	PAL	PAL	PAL
Audio Signal	Analog	Analog	NO	NO
Sensor Size	1/3″	1/3″	1/3″	1/3″
Sensor technology	CCD	CCD	CCD	CCD
Image format	4/3	4/3	4/3	4/3
Full Image size	-	-	-	-
Frame rate full resolution	50 frames/sec (analog mode)	50 frames/sec (analog mode)	50 frames/sec (analog mode)	50 frames/sec (analog mode)
Max Exposure time	-	-	-	-
ON board Memory	NO	NO	NO	NO
External Memory Card	NO	NO	NO	NO
External camera power	15V DC power supply	15V DC power supply	12V DC power supply	12V DC power supply
White Balance	Auto	Auto	Auto	Auto
Gain Control	Auto	Auto	Auto	Auto
Back light control	Auto	Auto	Auto	Auto
Exposure control	Auto	Auto	Auto	Auto
C-Mount connection	YES	YES	YES	YES
CS-Mount connection	NO	NO	NO	NO
Arm length	50 cm	65 cm	50 cm	65 cm
8mm objective	YES	YES	YES	YES

Accessories included:

All models (except VC-05 and C-E2): C-mount projection lens, 30 mm ring adapter, SCART for TV plug.

VC-05: C-mount projection lens, 30 mm / 30.5 mm ring adapters, SCART for TV plug. C-E2: 30 mm / 30.5 mm ring adapters, 1.8 m USB cable.

EYEPIECE CAMERAS - Educam & VC Series

C-E2



User-Friendly Eyepiece Microscope Camera

- Direct connection into the eyepiece tube (23 mm, 30 mm & 30.5 mm diameter)
- No additional adapters required
- Very useful for educational purposes
- Removable miniUSB cable
- Software included (for Windows): Optika Vision Lite, Optika LiteView and Optika ProView.
- Optika LiteView is also compatible with IOS & Linux.



Ready to use, replacing one eyepiece

USB / 4083.4	MIC / 4083.5	VC-05	C-E2
0.3 MP	NO	NO	2 MP
PAL 582 x 420	PAL 582 x 420	PAL 582 x 420	NO
PAL , USB2.0	PAL	PAL	USB 2.0
Analog	NO	NO	NO
1\3"	1\3″	1\3″	1\3.2"
CCD	CCD	CCD	CMOS
4\3	4\3	4\3	4\3
640 x 480	-	-	1600 x 1200
50 frames\sec (analog mode), 25 frames\sec (digital mode)	50 frames\sec (analog mode)	50 frames\sec (analog mode)	5 fps (1600x1200) / 7,5 fps (1280x1024) / 22fps (640x480)
Auto	-	-	Auto
NO	NO	NO	NO
NO	NO	NO	NO
15V DC power supply	12V DC power supply	12V DC power supply	PC USB
Auto	Auto	Auto	Auto
Auto	Auto	Auto	Auto
Auto	Auto	Auto	Auto
Auto	Auto	Auto	Auto
YES	YES	NO	NO
NO	NO	NO	NO
65cm	-	-	-
YES	YES	NO	NO

OPTICAL ADAPTERS - Adapters for OPTIKA Cameras and Microscopes

Upright/Inverted Microscopes (C-mount projection lens)



4

Cameras & Digital

(0.5x)

Upright/Inverted Microscopes (Focusable C-mount adapter)





M-620.1 (0.5x)

M-620.2 (0.65x)



Stereomicroscopes - SZM / SZO / SZP Series (Focusable C-mount adapter)



ST-090.1

(0.5x)

ST-090

(0.35x)



(0.65x)



ST-090.2

M-620.3 (1x)

		Upright			Inverted		Stereo		
		Monocular Binocular (Ø 23mm)	Trinocular (Ø 23mm)	Binocular (Ø 30mm)	Trinocular (Ø 30mm)	Trinocular (Ø 30mm)	Binocular (Ø 30.5mm)	Binocular (Ø 30mm)	Trinocular (Ø 30mm)
Camera models	Sensor size	Ecovision / B-60 B-150 / B-190 B-290 / B-380	B-190 / B-290	B-810 / B-1000	B-380 / B-510 B-810 / B-1000	IM-3 / IM-5	SFX	SLX / SZM / SZO SZP	SLX / SZM / SZO SZI
C-B1	1/3"	included with the camera	included with the camera	included with the camera	M-620	M-620	included with the camera	included with the camera	ST-090
C-B3	1/2"	included with the camera	included with the camera	included with the camera	M-620.1	M-620.1	included with the camera	included with the camera	ST-090.1
C-B5	1/2.5"	included with the camera	included with the camera	included with the camera	M-620.1	M-620.1	included with the camera	included with the camera	ST-090.1
C-B10	1/2.3"	included with the camera	included with the camera	included with the camera	M-620.1	M-620.1	included with the camera	included with the camera	ST-090.1
C-HA	1/2.8"	-	M-115	-	M-620	M-620	-	-	ST-090
C-HE	1/2.8"	M-115	M-115	M-115 + M-113.1	M-620	M-620	M-115 + M-113.2	M-115 + M-113.1	ST-090
C-HESC	1/2.8"	-	M-115	-	M-620	M-620	-	-	ST-090
C-HP	1/1.9"	M-114	M-114	M-114 + M-113.1	M-620.1	M-620.1	M-114 + M-113.2	M-114 + M-113.1	ST-090.1
C-HPSC	1/1.9"	-	M-114	-	M-620.1	M-620.1	-	-	ST-090.1
C-P3	1/2.8"	M-115	M-115	M-115 + M-113.1	M-620	M-620	M-115 + M-113.2	M-115 + M-113.1	ST-090
C-P6	1/1.8"	M-114	M-114	M-114 + M-113.1	M-620.1	M-620.1	M-114 + M-113.2	M-114 + M-113.1	ST-090.1
C-P6FL	1"	-	-	-	M-620.3	M-620.3	-	-	M-620.3
C-P8	1/2.5"	M-114	M-114	M-114 + M-113.1	M-620.1	M-620.1	M-114 + M-113.2	M-114 + M-113.1	ST-090.1
C-P20CC	1"	-	-	-	M-620.3	M-620.3	-	-	M-620.3
C-P20CM	1"	-	-	-	M-620.3	M-620.3	-	-	M-620.3
TB-3W	1/2"	-	included with the camera	-	M-620.1	M-620.1	-	-	ST-090.1
TB-5W	1/2.5"	-	included with the camera	-	M-620.1	M-620.1	-	-	ST-090.1
4083.F33	1/2.7"	M-115	M-115	M-115 + M-113.1	M-620	M-620	M-115 + M-113.2	M-115 + M-113.1	ST-090
4083.WiFi	1/2.5"	included with the camera	included with the camera	M-114 + M-113.1	M-620.1	M-620.1	included with the camera	included with the camera	ST-090.1

v 1.4.0 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice

Headquarters and Manufacturing Facilities

OPTIKA[®] S.r.I.

Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

Optika Sales branches

OPTIKA[®] Spain **OPTIKA**[®] China **OPTIKA**[®] India

spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com

OPTIKA[®] USA **OPTIKA**[®] Central America

usa@optikamicroscopes.com camerica@optikamicroscopes.com

OPTIKA SOFTWARE



OPTIKA SOFTWARE SUITES

OPTIKA SOFTWARE - Comparison chart

Software

• Before proceeding with the SW installation, please check the table below "Software Function list" to identify the most suitable software.



SOFTWARE FUNCTION LIST

	FUNCTIC	N	OPTIKA PRO VIEW	OPTIKA LITE VIEW	OPTIKA VISION LITE
	Simultaneous management of several cameras		Х	Х	Х
	GUI (Graphical User Interface)		Х		
	Report generator		Х		Х
	Archiving		Х	Х	Х
		Catalan	Х	Х	
		Chinese (simpl.)	Х	Х	
		Chinese (trad.)	Х	Х	
		Korean	Х	Х	
GENERAL	Language	English	Х	Х	Х
		French	Х	Х	Х
		German	Х	Х	Х
		Indonesian	Х	Х	
		Italian	Х	Х	Х
		Japanese	Х	Х	
		Polish	Х	Х	Х
		Russian	Х	Х	
		Spanish	Х	Х	Х
		Swedish			Х
		Thai	Х	Х	
		Turkish	Х	Х	

	FUNCTIO	N	OPTIKA PRO VIEW	OPTIKA LITE VIEW	OPTIKA VISION LITE
	Measurements on "live"		Х		
	Measurements on "captured"		Х		Х
		Line	Х		Х
S		Angle	Х		
Ë		Parallel lines	Х		
	2D Measurements	Rectangle	Х		
REMENT		Ellipse	Х		
L K		Circle	Х		
ASI		Annulus	Х		
MEASU		Arc	Х		
2		Curve	Х		
		Polygon	Х		
	Particle count		Х		
	Export to Excel		Х		Х

OPTIKA SOFTWARE - Comparison chart

SOFTWARE FUNCTION LIST

FUNCTI	ON	OPTIKA PRO VIEW	OPTIKA LITE VIEW	OPTIKA VISION LITE
Simultaneous management of seve	ral cameras	Х	Х	
IMAGE acquisition		Х	Х	Х
	tiff	Х	Х	Х
	jpg	Х	Х	Х
	bmp	Х	Х	Х
Image formats	png	Х	Х	
	рсх	Х	Х	
	jp2	Х	Х	
	dcm	Х	Х	
IMAGE acquisition		Х	Х	Х
	avi	Х	Х	Х
	wmv	Х	Х	X
	mp4	Х	Х	Х
VIDEO formats	asf	Х	Х	Х
VIDEO IOMALS	3gp	Х	Х	Х
	mov	Х	Х	Х
	h264	Х	Х	Х
	h265	Х	Х	X
Continuous automatic exposure		Х	Х	Х
Manual Exposure		Х	Х	Х
Mobile spot for exposure		Х	Х	Х
Resizable spot for exposure		Х	Х	Х
Colour acquisition		Х	Х	Х
Grey-scale acquisition		Х	Х	Х
Manual Time-Lapse		Х		Х
Automatic Time-Lapse		Х		
Fast Image Acquisition		Х	Х	Х
Focus Indicator		Х		
White Balance		Х	Х	Х
Black balance Background correction		Х		
		Х		
Dark Field Correction		X	Х	
Image Enhancement		Х	Х	Х
Live Histogram		Х	Х	Х
Flip	Horizontal	Х	Х	Х
	Vertical	Х	Х	Х
Rotate		Х		

Cameras & Digital

Shift Correction Х HDR (High Dynamic Range) Layer Management Х Х Text Overlay Х Ruler Overlay Х Measurement Overlay Х Х Х Grids

OPTIKA PRO VIEW

Х Х

Х

Х

Х

OPTIKA LITE VIEW

4083.Wifi, 4083.4 and 4083.EC2 work with Vision Lite only. Cameras with HDMI connection only, do not require any software.

FUNCTION

Several function of image processing (filters)

Colour Combine (Multi-Fluorescence Imaging)

Multiple image combining

PROCESSING

EDF (Extended Depth of Focus)

OPTIKA VISION LITE

OPTIKA Vision Lite - Extremely Intuitive Software

Optika Vision Lite has been designed and developed to be incredibly intuitive, simple and easy to use for customers needing a convenient solution to be combined with OPTIKAM cameras.

- » Friendly interface, multilanguage
- » Capture still images & stream live videos
- » Perform linear measurements
- » Export comprehensive reports

Friendly interface, multilanguage

Engineered for easy user interaction and optimized image acquisition, the main purpose of OPTIKA Vision Lite is ensure clear communication.

- •An efficient means to efficiently completing your jobs
- •Pleasant, easy-to-navigate menus
- Eight languages pre-installed, others upgreadable

Capture still images & stream live videos

Use the live preview to accurately focus your image and change parameters to obtain the perfect final result you are looking for. Images can be saved in different formats and even as test reports, including personal comments.

Additional features:

- Image stack acquisition
- Grid addition for rapid considerations
- Image flipping option available


Perform linear measurements


Perform linear measurements in an extremely way just by drawing a line after creating your preferred calibration based on the magnification.


- Accurate measurements through simple calibration
- Comprehensive data export (notes & measures included)
- Indicate particular objects in the image to add persona comments

Export comprehensive reports

Detailed test reports can be generated, printed and saved. Reports can be also customized with company logos.

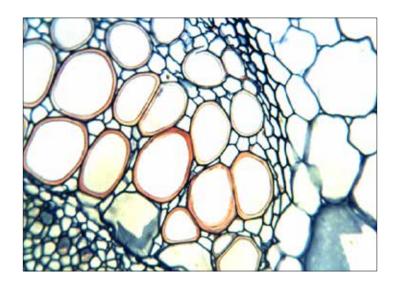

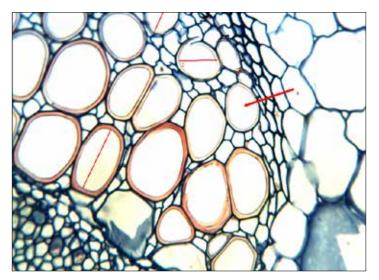

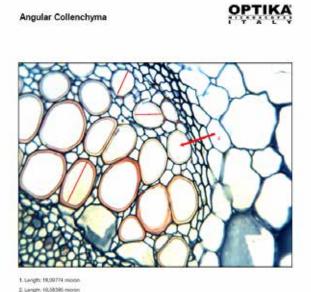

2. Length: 10.58395 micros 3. Length: 10.20079 micros 4. Indication

OPTIKA LITEView - Life is Easier

OPTIKA LITEView is a basic image acquisition software. The user who simply wants acquire a still image or a video, with no no need to perform measurements, has, with this powerful and intuitive software, the perfect solution.

- -) Simple management of «live» image
- -) Acquisition of still images or video
- -) Basic imaging functions
- -) Background correction



Cameras & Digital

Simple management of «live» image


Image preview is freely customizable by the user. A simple White Balance function with a mobile spot allows to perform the balance even on very small areas, once the specimen has been framed and focused.

Basic functions:

- Automatic or manual acquisition
- Possibility to have «live» and «capture» at different resolutions
- White Balance with mobile spot
- Background correction for the acquisition of perfectly illuminated images.

Capturing still images or video

Just select the option and the software performs: acquiring still images or videos is simply and intuitive.

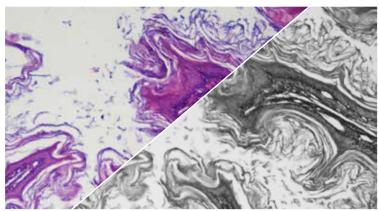

Record

Color / Grey scales

Basic imaging functions

Image parameters can be modified according user's needs. Color, Contrast and Gamma can be chaned in real time. More, it is possible to use a color camera in «SGrey Scales» modo in order to increase the camera sensitivity.

Background Correction

Any inhomogeneity of illumination of the microscope can be corrected by using the background correction function. This allows to obtain a faithful reproduction of the image without annoying inhomogeneity due to a not perfect illumination.



Camera List

Capture & Resolution

Snap

2592 x 1944

2592 x 1944

C-B5

Live:

Snap

No Background correction



With Background correction

OPTIKA PROView - Professional Image Analysis

OPTIKA PROView is a professional image analysis software. The user who needs to acquire an image or video and to perform a series of processings or measurements, can easily achieve incredible results thanks to this software. PROView incorporates all the functions of the LITEView package, but in addition allows:

- White Balance and Black Balance
- · Simultaneous management of several cameras
- Graphical User Interface fully customizable
- · Imaging of Multichannel Fluorescence Images with «pixel shift» function
- Multilanguage Software

Beginners? Experts?


An «On-line» manual will help any user (no matter on how expert he can be) to get the best from the software


Images always perfect

The management of the acquisition parameters allows to get always the best from your camera. White balance, black balance, background correction, «live» management of Colors, Contrast, Gamma, Gain and Exposure Time ensure to obtain a faithful image. A numerical focus indicator will ensure an optimal focusing, also on specimaens with different focal planes.

White Balance and Black Balance

It is possible to obtain the balance either on the whole frame or on a small ROI (Region Of Interest) of the image simply resizing and moving the spot in one part of the specimen

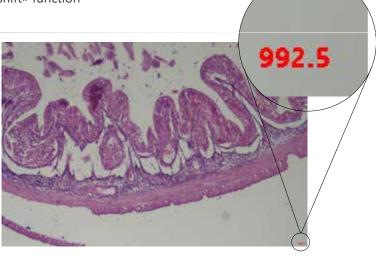

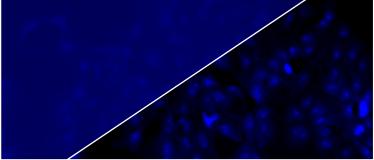


No black correction / Black correction

Multichannel Fluorescence Image processing


Acquire fluorescence images with a specific filtercube, use a false color for the used fluorochrome, get a single multichannel image is simply and intuitive.




DAPI (UV)

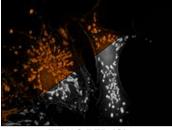

TEXAS RED (G)

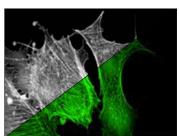

FITC (B)

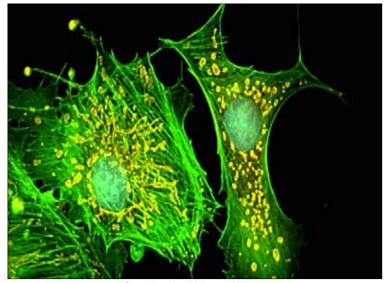

Combined multichannel image


OPTIKA PROView - Professional Image Analysis

«Pixel Shift» function

Fluorescence ilter cubes, sometimes, are not perfectly aligned.

During acquisition of multichannel luorescence images, this can cause a non perfect overlapping of the different signals, making the colocalization calculation almost impossible.

«Pixel Shift» function allows to correct these small misalignments:

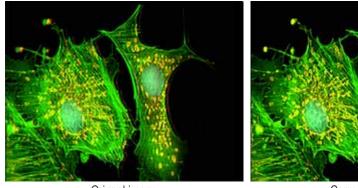


Orignal image

Corrected image

HDR (High Dynamic Range) acquisition

Acquisition of different images with different exposure times allows this function to create a final image where bright and dark zones of the specimen are perfectly displayed.

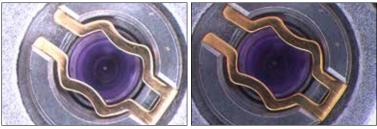

Standard Dynamic Range

High Dynamic Range

Extended Depth of Focus (EDF)

Acquire images with different focal planes, specially on specimens observed under a stereomicroscope, and to obtain a focused final image with a theoretical infinite focus. **EDF** function (also known as «Z-stack») allows a very refined image processing.




Single Focal Plane Images

EDF Image


Stitching & Tiling

Get an image with high resolution but, at the same time, have a wide view of the specimen under observation. Impossible? No. The multiple image alignment function allows to get a singe image starting from adjacent images of the specimen.

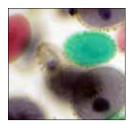



Separate Images

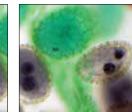


Stitched image


(4)

OPTIKA PROView - Professional Image Analysis

Measurements

User can perform measurements on the «live» image (no need to capture an image) and on captured images.

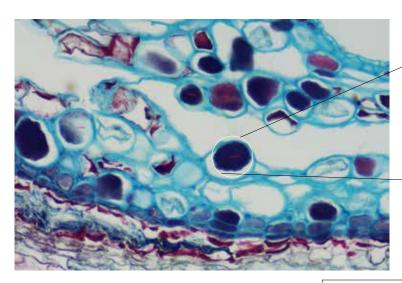

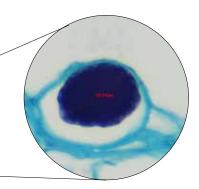

From Beginners To Experts


Measurements available:

- linear measurements
- angles
- circles
- annuli
- poligons
- touch count

Report Generator

At the end of the analysis it is possible to export images and measurement results either on a Excel sheet and on a Report Generator in MS Word format.

The template is freely configurable and can be modified according to laboratory standards.

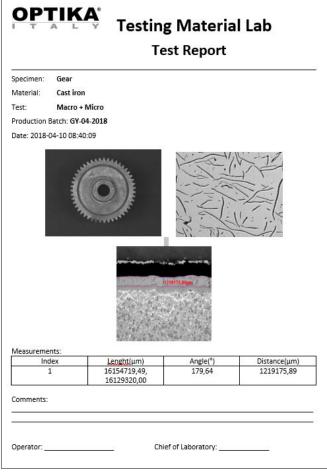

v 1.4.0 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

Headquarters and Manufacturing Facilities

OPTIKA[®] S.r.I.

Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

Optika Sales branches

OPTIKA[®] Spain OPTIKA[®] China OPTIKA[®] India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com

OPTIKA[®] USA **OPTIKA[®]** Central America

usa@optikamicroscopes.com camerica@optikamicroscopes.com

OPTISCAN




Cameras & Digital

OPTISCAN10 - 4083.SC10

CONVERT YOUR GLASS SLIDES INTO DIGITAL DATA !

Rapid and high resolution scanner to convert your slides into digital slides. The digital slide can be easily manipulated to see any location

- at any magnifications. Digitizing slides opens up a variety of new possibilities, like:
- Creating a database to be incorporated into a laboratory information system
- Networking slide libraries to be consulted from distant facilities and research institutes
- Sharing expertise for evaluation processes and discussing
- Information storing (digital data does not deteriorate, are secure from damages and losses)
- Main application fields are quality control & research, education, veterinary, histology / pathology, entomology / insectology, etc.

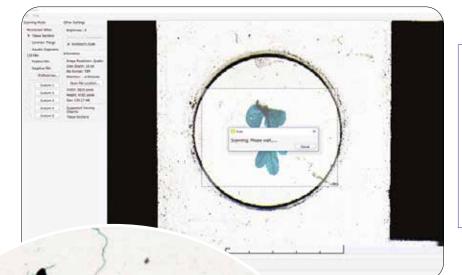

Main Features:

- High Resolution (up to 10.000 dpi)
- True & Neutral Color Fidelity
- White Balance & Distortion-free Images
- Dedicated Illumination (LED Transmitted Light)
- Efficient Scanning Area, Wide Field of View
- Impressive Scanning Speed (from 40 sec. to few minutes)
- High Sensitivity CCD Sensor
- Largest Field Of View, Better Than Any Camera

Ideal for:

- building up a comprehensive database of images for routine operations
- sharing expertise for evaluation processes
- archiving confidential patient information

OPTISCAN10 - Technical Specifications

OPTISCAN10 is an extremely convenient scanner for professionals, labs & teaching purposes, offering unmatchable price/performance ratio and coming along with a comprehensive but user-friendly software.

A ultra efficient, compact scanning device carrying high resolution features for spot detection with easy operation figure. It is equipped with a dedicated LED transmitted light system and high resolution CCD sensor, ensuring high sensitivity with low background noise.

Signal output	USB 2.0
Illumination	LED
Resolution	5'000 dpi (Normal), 10'000 dpi (Quality)
Allowed slide	Standard 24 x 75 mm
Scan view size	Any size, Max 24 x 36mm
Prescan function time	25 seconds
Scanning time (Normal)	1min 30sec (24 x 36mm); 40 sec (standard 15x15mm cover slide)
Scanning time (Quality)	2min 10sec (24 x 36mm); 1min (standard 15x15mm cover slide)
Always included	1.5 m USB cable, power supply, CD rom
System requirements	Windows XP service pack 2, Vista / win7 / win8 / win10 / 32-64 bit / USB 2.0
Supplied software	Multilanguage software for image scan
Capture features	Prescan, slide scan 24x36mm, crop scan, brightness, contrast, saturation, image flip



4

v 1.4.0 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

Headquarters and Manufacturing Facilities

OPTIKA' S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

Optika Sales branches

OPTIKA[®] Spain OPTIKA[®] China OPTIKA[®] India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com

OPTIKA[®] USA **OPTIKA**[®] Central America

usa@optikamicroscopes.com camerica@optikamicroscopes.com